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Massagepistol
Massage gun
Massaažipüstol
Массажный пистолет

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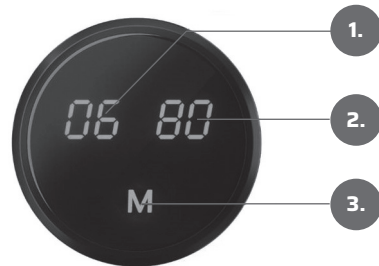
Danger

To reduce the risk of fire and personal injury, this product must be used according to the instructions below.

- For adults only.
- The device can only be used on dry and clean skin. The device should not be used on top of clothing.
- Gently press and move the device on the skin. Massage each area for about 60 seconds.
- Use only the vibrating head that gives the best results.
- If the pressure is too high, the massage can produce bruising. By signs of pain or discomfort the use should be discontinued immediately.
- Keep your fingers, hair or other parts of the body away from the shaft and back of the vibrating head, as hair may get stuck and skin may get pinched.
- Keep the vents of the device free from any objects.
- Do not immerse in water or allow water to enter the vents of the device.
- Do not drop the device or use it improperly.
- Charge only with a 24 V DC charger.
- Check the device and the used socket carefully before each use.
- Do not modify the device in any way.
- Never let the device run or charge unattended.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- The appliance is only to be used with the power supply unit provided with the appliance.
- The device contains a fixed battery (not replaceable).
- Children should be supervised to ensure that they do not play with the appliance.
- Children are not to clean or service the device without supervision.
- **WARNING:** Only use the charger provided with the device to charge the battery.

Use

Press the power switch on the top of the machine to turn the device on. The light turns on. Press the power switch again to turn off the device. The light turns off.

Press the power control button in the control panel to adjust the speed level. Press it once to set the vibrating speed on level one. Press the power control button twice to adjust the device to speed level two. A total of six levels can be used. Press the power control button for three seconds to pause the device. The speed level starts from level one when you press the power control button again. The number on the left side of the LED panel indicates the speed level (from 00 to 06). The number on the right side indicates the battery level. The battery level is based on the used speed level and the pressure applied to the muscle. That is why the battery level can go up and down during use.

Apply different levels of pressure for different muscles. When using the device for the first time, start with a low pressure and increase the pressure slowly. Use the device only on

muscles and other soft tissue on the body without causing pain or discomfort. Do not use the device on the head or any hard or boned part of the body.

If you have an injury or have had a surgery, always consult your doctor about using the device for recovery.

Replaceable vibrating heads

Change the vibrating head by pulling it out gently. The head should be easy to remove and using force should not be necessary.

We recommend trying different heads on different muscles to find a combination that suits you.

6 different heads and their use:

1. Large ball: big muscle groups
2. Small ball: small muscles
3. Fork: back muscles, calves, ankles, soles, and arm muscles.
4. Bullet: trigger points and palm muscles. Strong effect.
5. Flat: thighs, gluteal and pectoral muscles. For gentle massage and enhancing recovery.
6. Angular: abdominal muscles and larger muscle groups.

Charging

Charge the device for four hours before the first use. The device must not be used during charging.

When charging is in progress, the red LED light on the power button is lit. The light turns green when the battery is fully charged.

It is not recommended to completely drain the battery.

Make sure the power is off when you are not using the device. The operating time of the device depends on the speed level during use and the applied pressure.

WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance (XH1260-1000LG). The device has a built-in battery which must not be removed.

Warranty

- This product is intended for private use only.
- If there are problems with the product, contact the place of purchase. Present the receipt.
- The product warranty is valid for 12 months forward from the day the product is bought.
- The warranty does not cover normal wear and tear.
- The warranty does not cover damage caused by careless use or storage or acting contrary to the instructions.
- The warranty does not cover damage caused by an external factor (e.g. thunder or floods).
- The warranty ends if the product is modified in any way.
- During the validity of the warranty, the primary aim is to repair the product or replace it with a new one.
- The warranty is valid in the EU area.

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Manufactured for Broman Logistics Oy